

flexyfoot

Putting the spring back in your step!



- More grip
- Shock-absorbing air suspension
- All terrain use

What is Flexyfoot?

Shock-absorbing, anti-slip Flexyfoot is designed to replace the old fashioned ferrule and revolutionise walking aids. Flexyfoot gives improved grip on all types of terrain and eases the aches and pains associated with constant stress and impact on joints. They are suitable for people of all ages and levels of fitness.

Flexyfoot has been exceptionally well received by users and the company has received many unsolicited testimonials from delighted customers:

“I have been on crutches since 2007 and have a lot of arm and shoulder pain. Within two weeks of using Flexyfoot, all the pain disappeared”

Owen Reynolds

“My mum can't stop raving about her Flexyfeet and says they are the best thing she has bought...”

Becky Grant

“I could not resist trying Flexyfoot in the snow and ice and, boy, are they good! I feel more comfortable and confident on my crutches”

Fred Stevens, Beijing Paralympian archer

“Very impressed with the shock absorber – takes the pressure off when the stick touches the pavement. I take my hat off to whoever designed this”

Rita Burton

“Flexyfoot has made a great difference to my balance and comfort. I would certainly recommend them to anyone who has to deal with crutches for more than a few days”

Dr AB




- More grip
- Shock-absorbing air suspension
- All terrain use

How does Flexyfoot work?

The bellows is screwed into the 'collar' for easy fitting and replacement but the special thread also enables the foot to rotate without coming undone. This enables the user to twist and turn whilst being supported by their aid and reduces stress and strain to their wrist and shoulder joints.

Integral seal prevents ingress of dirt and water.

Patented bellows is sealed so that air acts as a spring when the bellows is compressed which cushions shock and impact.

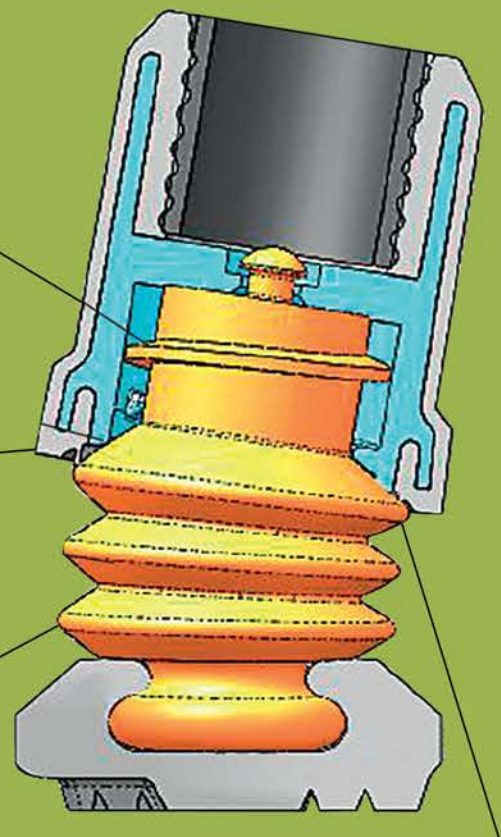

Conventional ferrule traps mud and stones. Flexyfoot's easy to clean tread helps to stay clear of debris making the transition from outside to inside easier and safer. It is also quieter!

The bellows allows the 'foot' to automatically adjust to the angle of the ground or floor, maximising the tread contact area. Grip is therefore greatly improved and wear is greatly reduced.

For replacement feet

- From your retailer
- Online at www.flexyfoot.com

flexyf^oo^t

Lancaster Healthcare Ltd t/a Flexyfoot

Unit 1, Brook Mill Industrial Estate,
Wrea Green, Preston, Lancashire, PR4 2PH

Collar guaranteed for manufacturing defects for 12 months.

540/Packaging/Hanger User Instructions 03-02. World-wide patents applied for.

